# Dr LUCAS

## NON-CONTACT THERMOMETER

**MANUAL  
ECT002**

# NON-CONTACT THERMOMETER


## IMPORTANT SAFETY INSTRUCTIONS

### READ CAREFULLY AND STORE THESE INSTRUCTIONS FOR FUTURE USE.

- This device has been designed for private use only and only for the purpose provided.
- Keep the device outside the reach of children.
- Make sure that the measuring lens is always clean.
- Put on the protective cap after each use.
- Do not immerse the device in water.
- Do not let the device fall down and avoid heavy impacts.
- Do not try to disassemble or modify the measuring device.
- Do not use the device if the measuring lens is damaged.
- If you hold the device in your hand for more than 3 minutes this may lead to wrong measuring results for any subsequent measurements. You should therefore avoid to hold the device in your hand for too long.
- If your child is younger than 2 years or if you do not measure your own temperature (e.g., the temperature of older persons) we recommend to take three measurements at the same distance. If the readings differ from each other choose the highest temperature.
- Minor deviations in repeated measurements are permissible.
- For the first use or when inserting new batteries wait about 10 minutes for the warm-up of the apparatus .
- Before measurement, make sure to remove hair and perspiration from the forehead.

## DESCRIPTION


- (1) LED display: It shows the measurement digits.
- (2) Backlight/laser key: It is in back light or laser light while turning off or turning on.
- (3) Setting key: It can be used for setting the device and the check the memories of the measurements.
- (4) ADJ key: It's used in checking memories of the measurements.
- (5) MODE key: It's used for measurement language selection between Chinese and English and the measurement target: surface or the body.
- (6) Laser light
- (7) Transducer
- (8) Measurement trigger: Pressing the trigger to measure.
- (9) Battery cover: Open to put or replace the batteries.
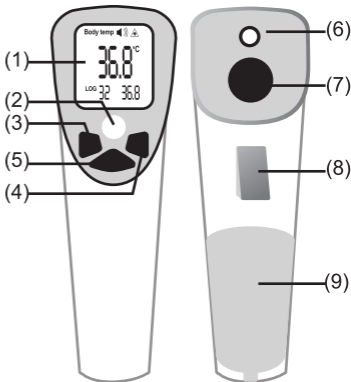
## PUTTING IN THE BATTERIES

The battery provided is a long-life alkaline battery (9V)

- Open the battery compartment, you will see the battery connector. See the “+/-” polarity as in the picture.
- Connect the battery “+” polarity to the “-” polarity on the battery connector. Push together to make sure that the connection is well done.
- Carefully put the connected battery **inside the** battery compartment. Then close the battery compartment.
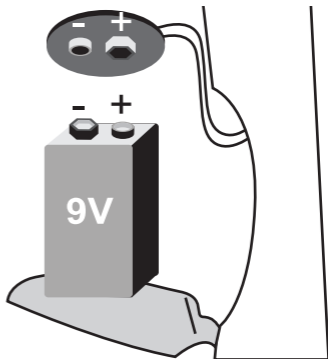
## BATTERY REPLACEMENT

It will show battery symbol on the left corner of the LCD display when the battery is low, please change the battery to make sure the device works properly.



## SETTINGS

• Press the measurement trigger and button MODE for 2 seconds, you can change the language between Chinese and English.


• press the measurement trigger and the button MODE again just once shortly, you can see the letter “surface” or “body” on the left screen of the LCD. Press MODE once to change between body and surface measurement.

• press the measurement trigger and then press SET key for 3 seconds till “F1” appears on the right corner of LCD, then press ADJ key to set the temperature unit between °C and °F.

• After set the temperature unit, press the SET key again, you will see F2 appears on the right corner of the LCD. F2 is for setting the temperature range, you can increase or decrease the temperature measurement range from -5°C to 5°C (±9°F) by pressing the ADJ key.

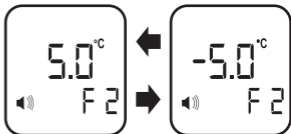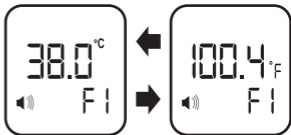- After set up the temperature range, press the SET key again, you will see F3 appears on the right corner of the LCD. You can select the sound setting by pressing the ADJ key "ON" means sound is on, and "OFF" means the sound is switched off.

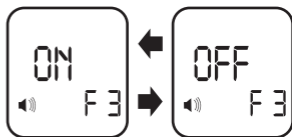- press the measurement trigger and the backlight button for 3 seconds to switch on or off LCD backlight.
- press the measurement trigger and the backlight button to switch on or off light, which will show you the point which you want to measure

## TEMPERATURE MEASUREMENT

Press the measurement trigger and straight the laser light at the surface or the body you want to measure. the lcd display shows the temperature.

### Surface temperature mode

The digits will continue updated if you are pulling the trigger to measure.

### Body temperature mode

Press the trigger to update measure the temperature and the current record will show on the lower right corner. Different temperature has different back light.

Green back light : below 37.5°C (99.5°F);

Orange back light: from 37.5°C to 38°C (99.5°F~100.4°F);

Red back light: above 38°C(100.4°F)

## OPERATIONS OF THE RECORDED MEMORIES


### Check the memories

Press the SET and ADJ keys to check the records

### Delete the memories

Press the ADJ key for 3 seconds to delete the current record and will show the previous record on the LCD display.(You will hear a sound "Du" when you switch on the sound function.)

Press the ADJ and SET keys at the same time to delete all the records.
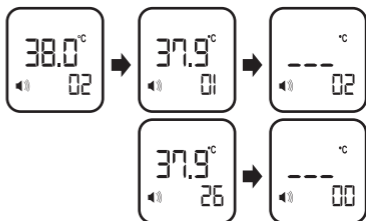
## SIGNS DISPLAY ON LCD

### SIGN "HI"

In the body temperature mode, when the temperature is over 43°C (109.4°) which is over the device measuring range. the LCD will display "HI".

### SIGN "LO"

In the body temperature mode, when the temperature is less than 32°C (90°F) which is over the device measuring range. the LCD will display "LO". This sign displays in various cases - please find below a list of the main cases

Reasons for "LO" display	Advice
Temperature measurement hampered by hair, perspiration.	Make sure that there is no obstruction prior to taking a temperature
Temperature measurement hampered by a strong air flux.	Make sure there is no air flux as this could interfere with the infrared system.
Temperature measuring are too often without pause, the device did not have the chance to boot itself.	Please make pause of at least 15 seconds between each measurement.
The measuring distance is too far	Please make sure the measuring distance is between 5 cm and 15 cm.

## MAINTENANCE AND STORAGE

- Use a dry cloth to clean the device after each use (recommended).
- Clean the measuring lens and sensor after each use.
- Do not use strong cleaning agents.
- If the measuring device is not used for a longer time, the battery should be removed. (Leaking batteries may damage the device).
- Protect the measuring device from direct sunlight, high temperatures, humidity and dust.
- Do not expose the device to extremely low (below  $-20^{\circ}\text{C}$  ( $-4.0^{\circ}\text{F}$ )) or high (above  $50^{\circ}\text{C}$  ( $122.0^{\circ}\text{F}$ )) temperatures.

## TECHNICAL PARAMETERS

Type of battery	9V
Storage temperature	$-10\sim 60^{\circ}\text{C}$ ( $14\sim 140^{\circ}\text{F}$ )
Usage temperature	$10\sim 40^{\circ}\text{C}$ ( $50\sim 104^{\circ}\text{F}$ )
Body Temperature measuring range	$32\sim 43^{\circ}\text{C}$ ( $90\sim 109^{\circ}\text{F}$ )
Body Temperature measuring accuracy	$\pm 0.3^{\circ}\text{C}$ ( $\pm 0.6^{\circ}\text{F}$ )
Surface Temperature measuring range	$0\sim 100^{\circ}\text{C}$ ( $32\sim 212^{\circ}\text{F}$ )
Surface temperature measuring accuracy	$\pm 1^{\circ}\text{C}$ ( $\pm 1.8^{\circ}\text{F}$ )
Minimum readings	$\pm 0.1^{\circ}\text{C}$ ( $\pm 0.1^{\circ}\text{F}$ )
Response time	0,5 s
Measuring distance	5~15 cm
Auto power off	7 s

## WARRANTY

1. This warranty is issued by Esperanza sp.j. Poterek with address:  
ul. Poznańska 129/133, 05-850 Ożarów Mazowiecki, Phone: 022-7213571 w.10 – called hereinafter Esperanza.
2. This warranty confers upon you specific rights and remedies against Esperanza. Please read this document carefully, by making a claim under this warranty, you will be deemed to have understood and accept its terms.
3. Esperanza warrants that any Esperanza product together with attached accessories is free from defects in materials and workmanship under conditions of its normal use.
4. If you make a claim under this warranty, Esperanza will, at its option, repair or replace any defective product into new one by using new or use factory-refurbished parts. In case of using factory-refurbished parts they will be equal to new products.
5. Defective product should be returned together with receipt, commercial invoice or any other proof of purchase to the place of purchase. It should include all accessories attached in original packaging. Esperanza has right to reject warranty if any accessories will be missed.
6. Esperanza is committed to respond to the complaint within 14 calendar days of its notification in the point of sale. Esperanza will repair the defective product in the shortest possible time. If will be need to import from abroad of any parts which are necessary to repair this time may be extended.
7. In case of replacing defective product onto new one, Esperanza has right to replace it with newer model with the same technical parameters if the old one model was replaced by new one during warranty period.
8. If you make claim under this warranty, you are deemed to have agreed that repair and (upon availability) replacement, as applicable, under the warranty services described herein, is your sole and exclusive remedy vis-a-vis Esperanza.

## EXCLUSIONS:

### This warranty does not extend to:

1. Any product not manufactured by Esperanza
2. Any product that has been damaged or rendered defective as a result of:
  - use it in other way than normal intended use
  - modification or internal interferes or changes of product
  - service or repair by anyone other than Esperanza service
  - improper store, transportation or packing
  - improper installation of product or software
  - mechanical, physical, chemical etc. damages
3. In no event shall Esperanza be liable for any damages, loss of business, profits, data or loss based on use Esperanza products with other devices.
4. Products with removed warranty and other labels and products spent or have hard use marks.

All products are covered by Esperanza's 2-year manufacturer's warranty while product is purchased by the consumer (the consumer is a person who buys product for the purpose of not professional or commercial activity). In the case of purchase for the company (for purposes relating to the activities of professional or commercial) warranty is 1 year from date of purchase. Esperanza products have all the necessary safety certificates required in the EU.

According to Directive 2002/96/EC, disposal of waste electrical and electronic equipment together with household waste is prohibited. Their components need to be recycled or disposed separately. Otherwise, they can cause the release of hazardous and noxious substances which pose a threat to the environment. The law requires each consumer to return unnecessary or waste electrical and electronic equipment to the producers, sales agents or waste treatment facilities. Detailed conditions are governed by national legislation. The above symbol placed on the product, in the user manual or on the packaging refers to this requirement. Thanks to proper segregation and compliance with the provisions relating to waste disposal each user contributes significantly to environmental protection.

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We, the manufacturer/importer (My, producent/importer):

**EDC Poterek Sp. Jawna**  
**ul.Poznanska 129/133**  
**05-850 Ozarow Mazowiecki**  
**Poland**

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Is in compliance with the essential requirements and other relevant provisions of the following directives:

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
**EMC 2004/108/C**  
**Year of issue (last 2 digits): 15**

The products is compatible with the following norms/standards:

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**EN 61326-1:2006**  
**EN 61326-2-2:2006**

Manufacturers/Authorised  
representative

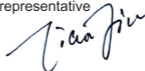Product Manager

Ozarów Mazowiecki 30-11-2015

CE certificate: [www.esperanza.pl/CE/ECT002.pdf](http://www.esperanza.pl/CE/ECT002.pdf)

Producent/ Producător/ Producer/ Výrobce/ Výrobca  
Productor/ Gamintojas/ Produtor/ Producteur/ Gyártó  
Produttore/ Produzent/ Παραγωγός/ Производитель:  
EDC Sp. Jawna. 05-850 Ożarów Mazowiecki  
ul. Poznańska 129/133 [www.esperanza.pl](http://www.esperanza.pl)  
MADE IN P.R.C.

© 2015 All rights reserved. Esperanza, Esperanza logo and Esperanza marks are owned by Esperanza and may be registered. All other trademarks are the property of their respective owners.