SMARTBMI SCALE BT

MT5513


User manual

EN

Introduction


Thank you for buying this modern scale with very interesting features. Device can work as simple scale, or smart scale when you use with smartphone – device will connect via Bluetooth and send data directly to smartphone (please download and install free application “OKOK” first). Application can show you BMI, body fat factor and many more. Please read this manual before using the device.

Overview

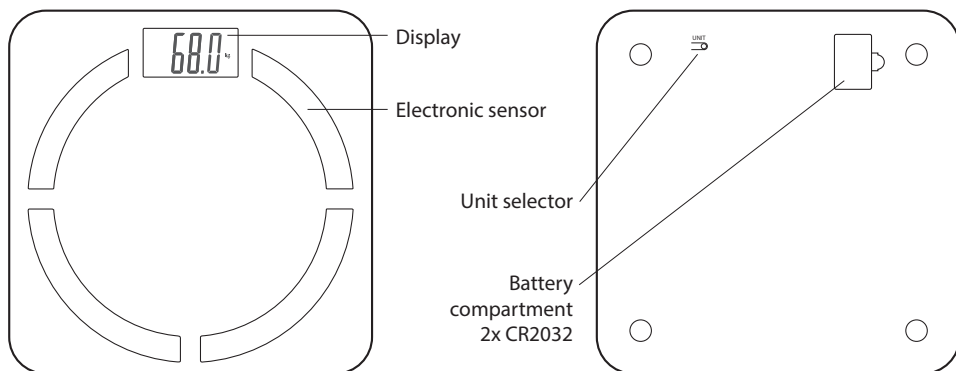Getting started

1. Install “OKOK” application (it is a free application available in GooglePlay™ for Android smartphone).

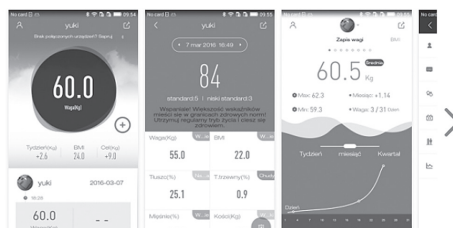2. Enable Bluetooth interface in your smartphone and stand on the scale with bare feet.
3. Scale will automatically connect with application.
4. Wait until for measurement is done.

Note! Scale lets a small amount of weak (less than 1mA) current flow through the human body, so as to detect the bioimpedance and estimate body fat. However, please be aware that anyone with an wearable or implantable medical electronic instrument, **such as a pacemaker**, must avoid using this device.

Due to continuous development specification and appearance of product are subject to change without prior notice. For technical support please visit www.media-tech.eu.

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