Oral-B®

powered
by **BRAUN**

PULSONIC

SLIM 2000-2900

English	6
Polski	9
Český	13
Slovenský	17
Magyar	21
Hrvatski	26
Slovenski	30
Lietuvių	34
Latviski	37
Eesti	40

UK

0 800 731 1792

IE

1 800 509 448

PL

801 127 286
801 1 BRAUN

CZ

800 11 33 22

SK

800 333 233

HU

(06-1) 451-1256

HR

01 / 66 90 330

SI

080 2822

LT

(8 5) 205 1272

LV

67798667

EE

667 5047

Internet:

www.oralb.com

www.braun.com

www.service.braun.com

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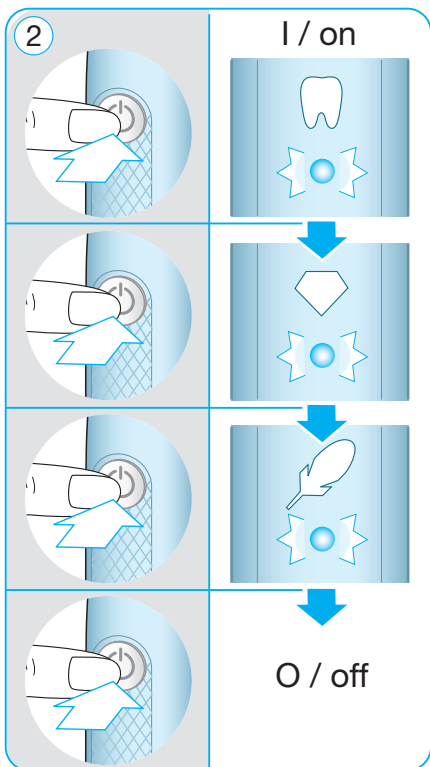
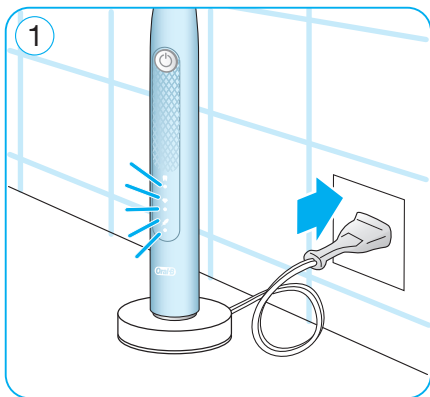
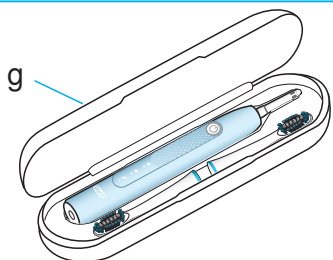
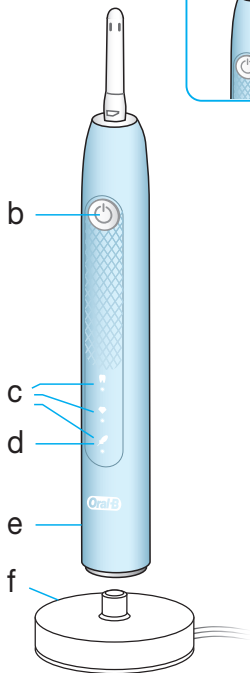
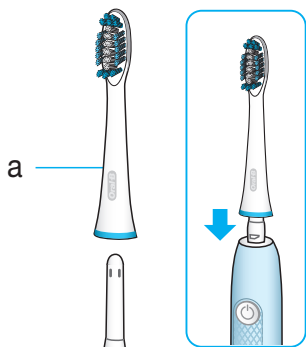
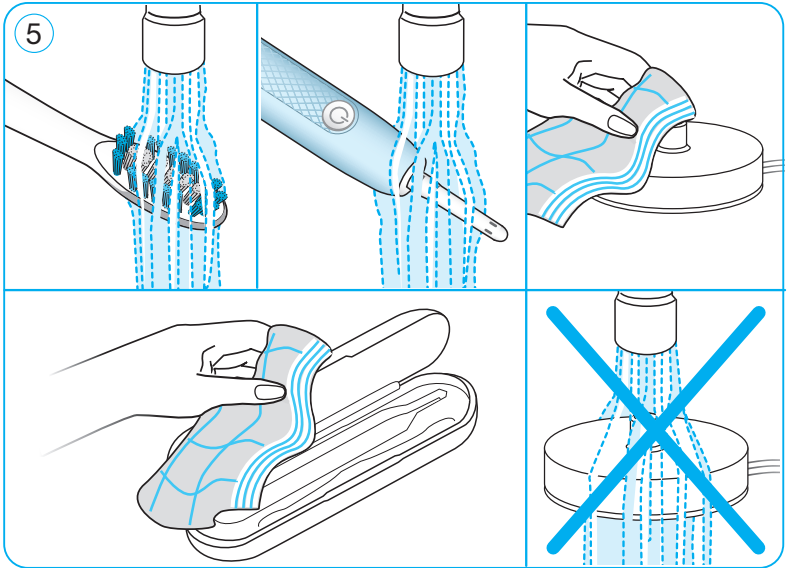
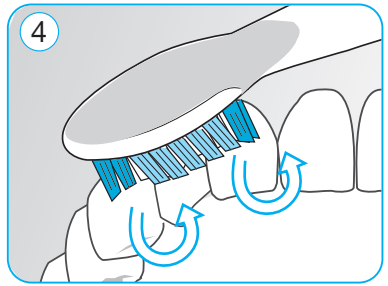
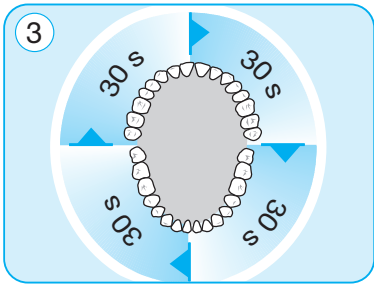
English

Welcome to Oral-B!

Your toothbrush has been carefully designed to offer you and your family a unique brushing experience that is both safe and effective. When using electrical products, however, basic safety precautions should always be followed.

Read all instructions before using.

WARNING

- Do not place the charger in water or liquid or store where it can fall or be pulled into a tub or sink. Do not reach for it when fallen into water. Unplug immediately.
- This appliance contains batteries that are non-replaceable.
- Do not disassemble the product except when disposing of the battery. When taking out the battery for disposal of the unit, use caution so as not to short the positive (+) and negative (-) terminals.
- When unplugging, always pull the power plug instead of the cord. Do not touch the power plug with wet hands. This can cause an electric shock.
- If you are undergoing treatment for any oral care condition, consult your dental professional prior to use.
- This toothbrush is a personal care device and is not intended for use on multiple patients in a dental practice or institution.
- It is recommended to use individual brush handles for each family member.
- This appliance contains magnets and generates magnetic fields. If you use a medical device or have an implanted medical device and you have questions about the safe usage of this toothbrush, consult the medical device manufacturer or a physician.

To avoid brush head breakage which may damage teeth:

- Prior to each use, ensure the brush head fits well. Discontinue use of the toothbrush, if brush head no longer fit properly. Never use without a brush head.
- If the toothbrush handle is dropped, the brush head should be replaced before the next use even if no damage is visible
- Replace brush head every 3 months or sooner if brush head becomes worn.
- Clean brush head properly after each use (see «Cleaning recommendations» para-

graph). Proper cleaning ensures the safe use and functional life of the toothbrush.

IMPORTANT

- Periodically check the entire product/cord for damage. Avoid sharp bending, jamming and squeezing of cord. If the product/cord is damaged, take it to a Braun/Oral-B Service Centre. A damaged or non-functioning unit should no longer be used. Do not modify or repair the product. This may cause fire, electric shock or injury. Never insert any object into any opening of the appliance.
- Usage by children under age 3 is not recommended. Toothbrushes can be used by children and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge, if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved.
- Cleaning and maintenance shall not be performed by children.
- Children shall not play with the appliance.
- Use this product only for its intended use as described in this manual. Do not use attachments which are not recommended by the manufacturer.

Only use charger provided with the appliance.

Description

- a Brush head
- b On/off button (mode selection)
- c Mode indicators
- d Charge indicator (with 2-mode device only)
- e Handle
- f Charging unit

Accessories (depending on model):

- g Travel case

Note: Content may vary based on the model purchased.

Specifications

For voltage specifications refer to the bottom of the charging unit.

Noise level: < 65 dB (A)

Charging and Operating

Your toothbrush is electrically safe and designed for use in the bathroom.

- Plug the charging unit (f) into an electrical outlet and place the handle (e) on the charging unit (picture 1).
- During charging the indicator light (d) fades in/out or all 3 mode indicators fade in/out (depending on model). Once the toothbrush is fully charged the light(s) turn(s) off. A full charge can take up to 24 hours and enables up to 2 weeks of regular brushing (twice a day, 2 minutes).
- If the battery is running low, the charge indicator (d) light is flashing or all 3 mode indicators are flashing for a few seconds after turning your toothbrush off (depending on model). It is recommended to recharge the toothbrush by then.
- To maintain the maximum capacity of the rechargeable battery, unplug the charging unit and fully discharge the handle with regular use at least every 6 months.

Using your toothbrush

Brushing technique

- Slide the brush head on the handle. Wet brush head and apply any kind of toothpaste. To avoid splashing, guide the brush head to your teeth before switching on your toothbrush.

- Place the brush head on the outside surface of your upper teeth. The toothbrush bristles should be placed against the teeth at a slight angle towards the gum-line.
- Turn your toothbrush on and start brushing in a slightly circular motion (picture 4). Apply light pressure during brushing. After a few seconds, guide the bristles to the next section. Brush inside, outside and chewing surfaces of your teeth with the same motion throughout the mouth for the recommended 2 minutes.
- To avoid splashing, turn your toothbrush off before removing from the mouth by pressing the on/off button (b).

You may experience a slight tickling or tingling sensation when you first use your toothbrush. This will subside as your teeth and gums will become accustomed to the cleaning experience.

In the first few days of using any electric toothbrush, your gums may bleed slightly. In general, bleeding should stop after a few days. Should it persist after 2 weeks, please consult your dentist or your dental hygienist.

Brushing modes (depending on model)

Your toothbrush offers different brushing modes for varying oral care needs:

- 🦷 **«Daily Clean»** – Gentle, effective mouth cleaning.
- 💎 **«White»** – Whitening mode for occasional or everyday use.
- 🦷 **«Sensitive»** – Gentle, yet thorough cleaning for sensitive areas

When pressing the on/off button (b), your toothbrush automatically starts in «Daily Clean» mode. To switch mode(s) successively press the on/off button. To turn your toothbrush off, press and hold the on/off button until the motor stops (picture 2).

Timer

A short stuttering sound at 30-second intervals (one break) reminds you to brush all four quadrants of your mouth equally (picture 3). A long stuttering sound (two breaks) indicates the end of the professionally recommended 2-minute brushing time. If the toothbrush is turned off during brushing, the elapsed brushing time will be memorized for 2 minutes. When pausing longer than 2 minutes the timer resets.

Pulsonic brush heads

The Pulsonic **CLEAN** brush head is recommended for everyday deep cleaning. Its bristles provide a thorough cleaning between your teeth, tooth surfaces and along the gumline.

The Pulsonic **SENSITIVE** brush head is recommended for gentle cleaning of sensitive areas. Its basic bristle trim and soft bristles clean thoroughly but gently.

Cleaning recommendations

Switch off the handle and remove the brush head. Rinse both brush head and handle thoroughly under running water. Wipe them dry before reassembling. Place the cleaned and dried handle on the charging unit. Charging unit (f) and travel case (g) should be cleaned with a damp cloth only. The charging unit should never be placed in water (picture 5). Store clean and dry tooth-brush/brush heads in the travel case only.

Subject to change without notice.

Environmental Notice

Product contains batteries and/or recyclable electric waste. For environment protection do not dispose in household waste, but for recycling take to electric waste collection points provided in your country

Guarantee

We grant 2 years limited guarantee on the product commencing on the date of purchase. Within the guarantee period we will eliminate, free of charge, any defects in the appliance resulting from faults in materials or workmanship, either by repairing or replacing the complete appliance as we may choose. Depending on availability replacement of a unit might result in a different color or an equivalent model.

This guarantee extends to every country where this appliance is supplied by Braun/Oral-B or its appointed distributor. This guarantee does not cover: damage due to improper use, normal wear or use, especially regarding the brush heads, as well as defects that have a negligible effect on the value or operation of the appliance.

The guarantee becomes void if repairs are undertaken by unauthorized persons and if original Braun/Oral-B parts are not used. To obtain service within the guarantee period, hand in or send the complete appliance with your sales receipt to an authorized Braun/Oral-B Service Centre.

Warranty replacement brush heads

Oral-B warranty will be voided if electric rechargeable handle damage is found to be attributed to the use of non-Oral-B replacement brush heads.

Oral-B does not recommend the use of non-Oral-B replacement brush heads.

- Oral-B has no control over the quality of non-Oral-B replacement brush heads. Therefore, we cannot ensure the cleaning performance of non-Oral-B replacement brush heads, as communicated with the electric rechargeable handle at time of initial purchase.
- Oral-B cannot ensure a good fit of non-Oral-B replacement brush heads.
- Oral-B cannot predict the long-term effect of non-Oral-B replacement brush heads on handle wear.

All Oral-B replacement brush heads carry the Oral-B logo and meet Oral-B's high quality standards. Oral-B does not sell replacement brush heads or handle parts under any other brand name.

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