
meriden®

Water Flosser Home & Travel

Before using this product, please read carefully this user's manual and keep it properly for later reference.

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
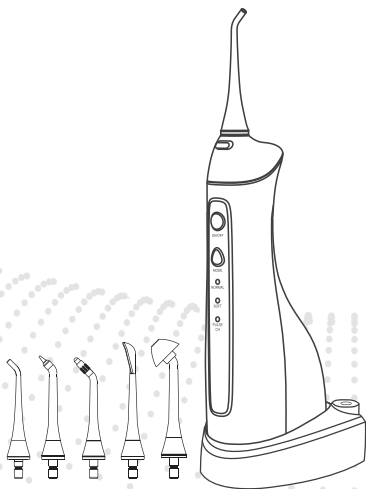
MD 1622

Safety Warning

To prevent burning, electrocution, fire or injury to people, take the following safety precautions:

- Do not plug this device into a voltage system different from the voltage system specified on the device or charger.
- Do not touch the charger with wet hands.
- Do not use the product if there is any damage to the cord, plug or the product itself, or the product is not working properly.
- Do not direct water under the tongue, into an ear, nose or other delicate areas, as the pressure may injure them seriously.
- Use this product only as indicated in these instructions or as recommended by your dentist.
- Fill the reservoir with water or dentist-recommended solutions only.
- Be careful not to drop the product. Never push any objects through the openings in this product. Power cord cannot be replaced- in the case of the power cord damage, the product must not be used.
- Keep the charging unit and cord away from heated surfaces.
- Do not use iodine, bleach or water insoluble concentrated essential oils in this product. The use of such substances may lead to impairment of the product performance and shortening of the product life.
- Remove any oral jewellery before using this product.
- Do not use this product if you have any open wounds on the tongue or in the mouth.
- If your physician advised you to receive antibiotic premedication before dental procedures, you should consult your dentist before using this instrument or any other oral hygiene aid.
- Children should be supervised to ensure that they do not play with the device.

Product description

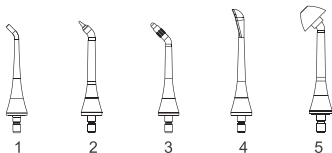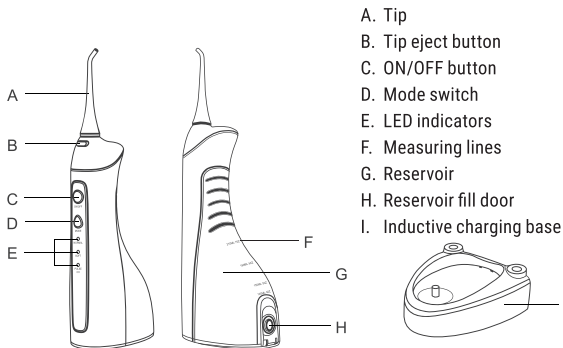
1. Classic Jet Tip - the tip removes debris from between the teeth. Intended for general everyday use.
2. Pocket Tip - thanks to its rubber finishing gently and thoroughly cleans periodontal pockets and the gum line.
3. Orthodontic Tip - thoroughly flushes areas between the teeth additionally cleaning areas of dental braces locks and archwires.
4. Tongue Cleaning Tip - for cleaning of the tongue - removes bacteria and food debris from the tongue.
5. Nasal Tip - flushing nasal passages due to the common cold, flu or allergies.

How to use


Charging the unit

1. Charging indicator light will turn red and blink when the battery is in the need of charging.
2. Plug the power cord into the wall socket (100-240V, 50/60Hz) and place the water flosser on the inductive charging base to charge the unit.
3. While charging, indicator on the handle will flash and when fully charged (approximately 24 hours) it will stay on continuously.
4. A full charge will last more than 30 minutes of standard use.

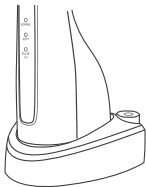Installing and removing the tip


- Insert the tip into the center of the knob at the top of the flosser handle.
- Press firmly until the tip clicks and locks into place.
- To remove the tip from the handle, press the tip eject button and pull the tip from the handle.

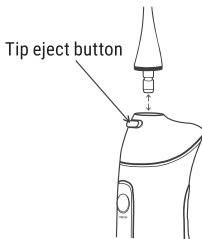Caution: Do not press eject button while the unit is running.


Filling the reservoir

1. Flat on the water flosser keeping the reservoir fill door upwards.
2. Open the reservoir fill door and fill the reservoir with lukewarm water.
3. Close the reservoir fill door when the reservoir has been filled.

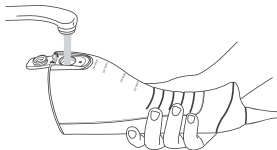Adjusting the pressure settings

- You can adjust the pressure setting on your water flosser by pressing the mode button on the handle.
- Start with the lowest setting (soft<normal<pulse) and increase to your liking.

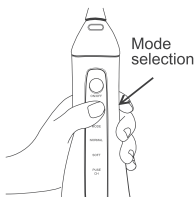
Recommended cleaning tips

- Lean over and place the tip in your mouth. Aim the tip at your gumline and turn the unit ON.
- For best results, start with the back teeth and work towards the front teeth.
- Glide the tip along the gumline and pause briefly between the teeth until you have cleaned the inside and outside of both the upper and lower teeth.
- Direct the jet stream at a 90° angle against your gumline. Slightly close the lips to avoid water splashing out but at the same time allow it to flow freely from your mouth into the bathroom sink.
- For best results, keep the unit in an upright position during use.

After use

1. Turn the unit off.
2. Empty any remaining liquid in the reservoir.
3. Clean the unit directly under the tap as it is IPX7 water-proof design.
4. Dry the unit with a cloth before putting it back on the charging base.

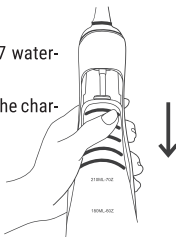Cleaning

- When necessary clean the product with use of a soft cloth and mild non-abrasive cleanser to wipe the product exterior.
- Water reservoir is removable for easy cleaning.
- It is recommended that that water reservoir is removed and cleaned once a week.

Storage

- Store the device out of children's reach.
- Storage space should be dry and clean with no direct sun exposure and away from heat sources.
- If the water flosser needs to be stored for an extended period of time, (over 3 months) be sure to fully charge it.

Troubleshooting

| Problem | Possible cause | Solution |
|---|--|---|
| Does not work | The product is new - it has just been purchased or it has been idle for more than 3 months | Charge the device for 24 hours. |
| It only works for a few minutes even after charging | Charging time was insufficient | Charge the unit for at least 24 hours |
| | The battery has reached the end of its life cycle. | Contact seller for support. |
| Water pressure is very low | Dirt has stuck in the tip. | Replace the tip |
| | The tip is deformed | |
| | Dirt has stuck in the filter | Clean the filter |
| Water does not come out | The reservoir is empty | Fill up reservoir |
| | The water irrigator has been tilted too much when used | Use the device in the vertical position |

Environmental Protection

If you want to discard of this product, please, do not place it together with domestic waste, but send it to the officially designated recycling centre and dispose of the battery at the officially designated recycling centre. Doing so can support environmental protection.

Warranty and Service

The product comes with 24-month warranty.

Explanation:

Please fill out the warranty card and keep it for future use. This product can enjoy two years of free warranty service counting from the date of purchase. Free warranty service does not apply in the following case of:

1. Damage caused by improper use or improper storage;
2. Self-disassembly or repair of the product outside of MERIDEN's special repair centers;
3. Lack of valid proof of purchase or the valid proof of purchase was altered;
4. Product damage due to voltage instability or voltage exceeding specified range;
5. Product exceeded warranty period or scope of warranty

After-sales service card:

Product model:

Date of purchase:

Customer name:

Contact information:

Fault description:

Warranty covers the unit and the charger (excluding the jet tip).

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