BEAUTIFLY

B-Lumix

P R O

series

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Specifications – B-Lumix PRO

Size	56x187x45.8 mm
Weight	Approx. 154 g
Rated power	10W
Rated voltage	DC 5V
Rated current	Approx. 2A
Timer	6 min
Working mode	RF (1000 KHz) / EMS 5-9µA
Battery	7.4V polymer lithium battery; 600mAh

Introduction – one step towards firm, radiant and rejuvenated skin


The loss of skin firmness is affected by a number of factors: UV radiation, genes, lifestyle, gravity, loss of fat, or hormones. The loss of collagen and elastin, the fibres of which form a scaffolding for the skin, is primarily responsible for the formation of wrinkles. Collagen gives your skin the right tone, elasticity and flexibility. The continuity of the skin's cell renewal processes also decreases with age. As a result, the skin will lose the proper level of moisture and good condition.

If your skin is deprived of vitality, it is giving in to the passage of time and showing more or less visible signs of aging, the B – Lumix PRO is the perfect answer to your problems. Skin rejuvenation with the use of radio waves (RF) and EMS currents is the best non-surgical facelift for your skin. Complete skin regeneration, dramatically improved vitality, firmness and youthful appearance in one!

How it works

The B-Lumix PRO maximises the production of collagen by means of thermal power combined with EMS micro-currents for intense facial skin rejuvenation.

1. 24-carat gold-plated energy head
2. Epidermis
3. Dermis
4. The B-Lumix PRO heat output is concentrated deep in the dermis
5. Heat passing into deeper layers of the skin
6. Subcutaneous tissue
7. Surface neuromuscular system

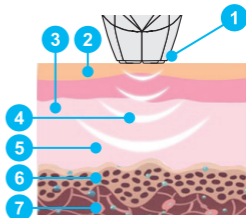Radio Frequency (RF)

Radio Frequency (RF) is a proven technology used in the medical and cosmetic industries. At a strictly controlled temperature, a magnetic pulse acts directly on the cells, producing thermal energy. As a result of the heating of collagen fibres, they shrink, resulting in an immediate effect on the skin. At the same time, the radio frequency activates the production of collagen. The results are visible shortly after the treatment: face oval lifting, skin firming, wrinkles shallower.

A series of treatments (2-3 times a week for 6-8 weeks) produce a long-term effect on skin structures, removes wrinkles, reduces "bags" under the eyes, restores youthful facial features.

Facial lifting massage using EMS (anti-aging) currents

This type of massage quickly firms the face and restores elasticity to the skin. It smooths wrinkles (skin-ironing effect), effectively improves the facial contour, and reduces expanded pores. EMS micro-currents stimulate skin cell activity, so intense regeneration occurs even at deep skin levels. EMS massage stimulates blood microcirculation, resulting in reduced discoloration, improved skin tone and good hydration. The skin is nourished, oxygenated and full of vitality.

24-carat gold-plated energy head


To achieve the best possible non-surgical lifting results, the instrument head is covered with 24-carat gold. Micronized gold increases collagen stimulation, making it easier to achieve lasting effects quickly. Pure gold is anti-bacterial – it reduces irritation and inflammation. The energy head is a state-of-the-art technology that makes the treatment effective, yet completely painless

Instrument specifications B-Lumix PRO

1. Power button
2. RF function indicator
3. EMS function indicator
4. Charging base

The B-Lumix PRO set consists of the following components:

5. Product
6. Satin bag
7. USB cable
8. Operator's manual

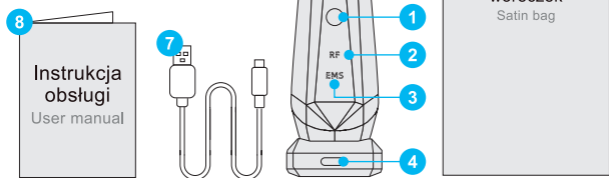Skin care – rules for carrying out treatment

To condition your skin with the B-Lumix PRO:

1. Wash your face.
2. Put an appropriate amount of gel/serum or other cosmetic product with slippery texture on the area to be used (face or neck).
3. Hold in the button on the instrument to start it.
4. Select the desired mode.

Important! During the first treatment session, it is recommended to set it at low power and then adjust the power level to the skin tolerance for subsequent treatments.

It is recommended that you perform a skin sensitivity test before using the instrument for the first time. The test must be performed 24 hours before the first use of the B-Lumix PRO. To do this, apply a small amount of gel/serum (or other cosmetic product with a slippery texture) to the inner part of the arm and treat the area with the B-Lumix PRO.

Make sure your skin does not react negatively to gel/serum (the cosmetic product you use) or the treatment.

5. Divide your face into 4 treatment areas:
 - Right cheek
 - Left cheek
 - Forehead
 - Outside corner of the eye
6. Place the instrument on your skin, make sure you lightly press the instrument in order to ensure that the energy head is fully in contact with your skin.
7. Start moving the instrument with very small circular movements on a small part of the treatment area (approximately twice as large as the treatment head). The energy of the instrument will be released into the skin, providing the skin comfort and warmth.
 Important note: Be sure to move the instrument slowly over the skin continuously to ensure that the treatment area receives adequate uniform energy. Do not hold the instrument in one position without movement.
8. Continue to work on the entire treatment surface slowly in a circular motion for 2-4 minutes. Maintain consistent energy by massaging the entire treatment area.
 Note: The effect of the procedure can be enhanced by performing "lifting" movements from the jaw line upwards toward the cheeks.
9. After 2-4 minutes of treatment, you can move to the next treatment area repeating steps 6-8. Switch to RF + EMS mode (refer to "Modes" in the instructions). It is recommended to start with a low power level.

Important: Use RF + EMS ONLY in the area of the jaw line. Do not use RF + EMS on other areas of the face.

10. Start moving the instrument by making very small linear movements along the jaw line.
11. For optimal results, use the instrument for 4 minutes on each side. After the whole treatment area has been massaged, the treatment session is ended. A natural effect of the treatment may be reddening of the skin (erythema) and a slight swelling in the treatment area. Both should disappear after a short time.
12. Hold in the power button to turn off the product.
13. If necessary, you can apply a moisturising cream.
14. Clean the energy head of the instrument with a damp cloth or paper towel, and place the instrument on the protective stand to protect the energy head for the next use.

Warning: Do not immerse the B-Lumix PRO in water.

Safety mechanism

The automatic temperature optimisation mechanism in the instrument protects your skin against overheating.

Treatment recommendations

The following treatment guidelines are recommended when using the B-Lumix PRO:

- It takes 8-12 minutes to treat the whole face. To reduce treatment time, focus only on the selected areas you want to treat. Treatment times vary per person.
- For optimal results, use the B-Lumix PRO 2-3 times a week for 6-8 weeks.
- A break of at least one day should be taken between treatments in the same treatment area.
- In order to maintain long lasting effects after the complete treatment, we recommend using B-Lumix PRO 2-4 times a month or as needed.

How to use the instrument

1. Charging

- 1) When you use the instrument for the first time, charge the instrument to make sure that the battery is fully charged (if the instrument is fully charged, there is no need to charge it).
- 2) When charging, place the instrument on the charging base and insert the USB connector into the charger plug and connect it to the general household power source. The LED indicator will light red during charging and the red light will go out when fully charged.
- 3) During the charging process, the 415 nm wavelength lamp disinfects the treatment heads for 5 minutes.
- 4) The charging time is approximately 120 minutes and the operating time is approximately 50 minutes.

2. Starting and stopping

Starting: press and hold to turn on the power in the off state. The sound indicator will emit the sound "D, D" when the power is turned on and off.

3. Modes

In start-up mode, briefly press to switch between modes from Mode 1, Mode 2, Mode 3, Mode 4, Mode 5. The sound indicator will emit "D" and the LED indicator will change colour.


Mode 1: RF Treatment (Level 1), the LED indicator will be continuous green.

Mode 2: RF Treatment (Level 2), the LED indicator will be continuous blue.

Mode 3: EMS procedure (Level 1), the LED indicator will be continuous green.

Mode 4: EMS procedure (Level 2), the LED indicator will be continuous blue.

Mode 5: RF and EMS treatment (the LED will light up with a continuous blue light during automatic switching between RF and EMS).

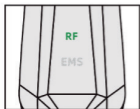
(Mode 1)


(Mode 2)

(Mode 3)

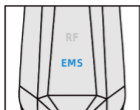(Mode 4)

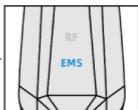(Mode 5)

Temperature sensor

The treatment heads are equipped with a precise NTC temperature sensor.

When an RF treatment is started, the sensor detects skin temperature; if the temperature is below 42°C, the RF treatment will begin to operate. If it is above 42°C, RF treatment will be interrupted.

Low battery level reminder

When the instrument battery is too low, the sound indicator will sound "D" every 5 seconds.

Automatic time measurement

The automatic time measurement function sets the time to 6 minutes. After this time, the product stops working and must be restarted to continue using it.

Safety

Please read the complete user manual carefully before using the product. Always observe all warnings and recommendations. When using B-Lumix PRO, all of the following precautions must be observed.

Warnings and recommendations

The B-Lumix PRO should be used correctly according to the instructions given in this manual. To use the instrument safely, observe all of the following warnings and precautions:

1. Do not place or store the instrument near water sources, such as a sink or bathtub.
2. Do not immerse the instrument in water.
3. Do not leave the instrument unattended when it is plugged in.
4. Before use, check the instrument and its adapter for visible damage. Do not use the instrument if the electrical cord is damaged or if the instrument is disassembled.
5. Only use the instrument with the electrical cord and adapter supplied.
6. Keep the instrument out of the reach of children.
7. Make sure your hands are completely dry when you operate the instrument.
8. Switch off the instrument after use or before cleaning and maintenance.
9. Use gel/serum (or any other cosmetic product with a slippery texture) to ensure proper lubrication during the procedure, good RF/EMS conductivity, and better moisture in the skin.

Safety rules

1. Use the instrument **only** for its intended purpose as described in this manual.
2. Do not use the instrument on **eyelids**, near the eye, in the thyroid area, on the **lips**, in the ear, on the breast or near the **genital** area.
3. The **feeling** of using the instrument should be nice and **pleasant**. If you experience discomfort, stop using the instrument immediately.

DO NOT USE THE INSTRUMENT IF ANY OF THE FOLLOWING CONDITIONS APPLY TO YOU:

1. You are under 18.
2. You have a cardiac pacemaker or an internal defibrillator.
3. You have **metal implants** in the treatment area, such as gold threads. This does not apply to dental fillings or implants, which do not present problems.
4. You **currently** have cancer or have had cancer in the past, especially skin cancer or pre-cancerous moles.
5. You have serious co-existing conditions, such as cardiac disorders.
6. You are pregnant or breastfeeding.
7. You have a compromised immune system due to immunosuppressive disease such as AIDS, or you are using immunosuppressive drugs.
8. You have any heat-stimulated **illness**, such as recurrent herpes simplex (in the treatment area), among others.
9. You have a history of keloids, poor wound healing or delicate skin. In the **last three months**, a surgical procedure, a laser procedure or a deep chemical peeling has been performed in the treatment area and you are **still** in the process of complete healing.
10. You have permanent **fillers** in the treated area.
11. If there are any temporary dermal fillers in the treatment area, **consult** your healthcare professional before using the instrument.
12. In the last three months, Botox has been administered in the treatment area.
13. You have suspicious changes in the treatment area.
14. You are **currently** using medications, herbal preparations, vitamins, or food supplements that may cause increased sensitivity or dry skin

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Electrical and electronic devices as well as batteries and accumulators must not be mixed with other household waste. This is indicated by the symbol of the crossed out basket on the packaging. Worn-out electric and electronic equipment as well as batteries and accumulators should be delivered to local separate waste collection points or to the seller, and for details inquire in your commune. Electrical and electronic equipment as well as batteries and accumulators may contain hazardous substances, mixtures and components harmful to the environment and human health. Improper handling or damage to them may lead to health detriment or environmental pollution during further disposal or recycling of the equipment.

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Beautifly Sp. z o.o., Plac Bankowy 2, 00-095 Warszawa

English

Hereby, Beautifly Sp. z o. o. declares that the device equipment Beautifly B-Lumix PRO is in compliance with Directive 2014/30/UE.

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TERMS OF WARRANTY

1. The product warranty is valid for a period of 24 months.
2. The warranty covers the period from the date of receipt of the device from Beautifly Sp. z o. o. and covers defects resulting from causes inherent in the received device.
3. The device will be repaired as soon as possible, not exceeding 21 working days from the date of acceptance of the device for repair by an Authorized Service Center.
4. The customer has the right to apply for replacement of equipment free from defects, if the Authorized Service finds that the removal of the defect is impossible or requires excessive costs.
5. If only part of the product is defective and can be detached from the part of the product working properly, the Customer's rights resulting from this warranty limit only to the defective part of the product. In order to benefit from the warranty conditions, please send a complaint notification together with the warranty card to the e-mail address serwis@beautifly.eu

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