# ORDO

## SONIC+

ELECTRIC TOOTHBRUSH

Important Information

Wichtige Informationen

Une information important

Información importante

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Önemli bilgi

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Viktig information

Vigtig information

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**ordolife.com**

Ordo, Manchester, WA14 4DR, UK.

Authorised Rep Compliance Ltd., Ground Floor,  
71 Lower Baggot Street, Dublin, D02 P593, Ireland.

# English

## Recycling

At Ordo we are always looking for solutions to reduce packaging and product materials that cannot be recycled effectively. We ask every member of the Ordo family to help us by recycling responsibly.

For more information, please visit [ordolife.com/recycling](http://ordolife.com/recycling)

## Troubleshooting

Problem	Possible Cause	Solution
The toothbrush isn't charging.	The charging base may be defective.	Try other sockets or USB ports and charge the toothbrush for at least 8 hours.
	The USB connector may be defective.	
	The charging base is defective.	Replace the charging base within the warranty period. Contact our customer service team.
The toothbrush operates with only a small amount of vibration.	The rechargeable battery is depleted.	Charge the toothbrush for at least 8 hours.
The toothbrush does not run for 2 minutes.	The rechargeable battery is depleted.	Charge the toothbrush for at least 8 hours and make sure that the toothbrush is properly placed on the charging base.

## Technical Data

Model number: SP2000-CG, SP2000-RG, SP2000-WS, SP2000-PV, SP2000-MG, SP2000-IB

Charging Base

Input: DC 5V, 1.5W

Protection class: III

IP code: IPX7

Dimensions (LxWxH): 6.2cm x 6.2cm x 4.3cm

Toothbrush

Rechargeable battery: 3.7V, 700mAh

Operating time: 56 cycles/full charge

Timer: 0.5 second pause after 30 seconds, deactivation after 2 minutes.

Dimensions (LxWxH): 2.6 x 2.7 x 25.4cm

## Warranty

This product is covered by a 2 year warranty. The extent of the warranty on this product is limited to technical defects caused by faulty production processes. In the event you want to claim warranty, please be sure to contact our customer service team for instructions. For more information please visit [ordolife.com/warranty](http://ordolife.com/warranty)

## 60 Day Money Back Guarantee

Satisfaction guaranteed – if you don't love your brush within 60 days we'll give you a refund. For more information on our guarantee please visit [ordolife.com/guarantee](http://ordolife.com/guarantee)

## Customer Service

For more information on Ordo's products and services please visit [ordolife.com](http://ordolife.com)

If the device is damaged, defective, in need of repair or you require assistance please contact our customer service team:

Email: [hello@ordolife.com](mailto:hello@ordolife.com)

Support: [help.ordolife.com](mailto:help.ordolife.com)

Our team are always happy to help.

## WARNINGS & CAUTIONS

### Safety Instructions

The toothbrush is designed exclusively for cleaning teeth. It is only intended for private use and is not suitable for commercial purposes. Only use the toothbrush as described in these operating instructions. Any other use is deemed improper and may result in damage to property or even personal injury. The toothbrush is not a child's toy. Ordoclub Ltd accepts no liability for damage caused by improper or incorrect use.

### Risk of electric shock

A faulty electrical installation or excessive mains voltage may result in an electric shock.

- Only connect the charging base if the mains voltage of the socket corresponds to the specification on the nameplate.
- The main cord of the charging base cannot be replaced. In the event of damage, the charging base must be disposed of.
- Only connect the charging base to an easily accessible socket so that you can quickly disconnect it from the mains in case of a failure.
- Do not operate the toothbrush and the charging base if they exhibit visible damage or the main cord or USB connector is defective.
- Do not open the housing; instead, have a qualified professional perform any repairs. Contact a qualified workshop for this. Liability and warranty claims are waived in the event of repairs performed by the user, improper connection of the device or incorrect operation.
- Only parts that comply with the original device information may be used for repairs. Electrical and mechanical parts, which are essential for providing protection against sources of danger, are located in this toothbrush.
- The charging base must not be operated with an external timer or with a separate tele-control system.
- Do not immerse the charging base or the main cord in water or other liquids.
- Never touch the USB connector with damp hands.
- Never pull the USB connector out of the socket by the main cord.
- Never use the main cord as a carrying strap.
- Keep the charging base, toothbrush and USB connector away from open flames and hot surfaces.
- Do not kink the main cord and do not lay it over sharp edges.
- Only use the charging base in rooms located indoors.
- Never setup the charging base such that there is a risk of it falling in a bath tub or sink.
- Never attempt to grab the charging base if it has fallen into water. In such a case, immediately disconnect the USB connector.
- Only charge the toothbrush with the enclosed charging base.
- Make sure that children do not insert any objects into the charging base or the hand piece.
- Disconnect the USB connector when you are not using the

charging base, after the charging process, before you clean the toothbrush or in the event of a fault.

### **Risk of chemical burns**

Fluid leaked from the rechargeable battery may lead to chemical burns.

- Do not allow acid from the rechargeable battery to come into contact with skin, eyes or mucous membranes.
- In the event of contact with acid from the rechargeable battery, immediately flush the affected areas with plenty of clean water and seek medical help.

### **Danger of explosion**

Improper handling of rechargeable batteries may cause an explosion.

- Never short circuit rechargeable batteries.
- Never dismantle rechargeable batteries.
- Never throw rechargeable batteries in fire.
- Never expose rechargeable batteries to excessive heat or direct sunlight.
- Do not in any case open the rechargeable batteries.

### **Danger for children, persons with impaired physical, sensory or mental capacities**

e.g. partially disabled persons, older persons with reduced physical and mental capacities, and those lacking experience and knowledge (e.g. older children).

- This toothbrush may be used by children aged 8 years and above, persons with impaired physical, sensory or mental capacities, and those lacking experience and knowledge if they are supervised or have been instructed in how to safely use the toothbrush and have understood the risks associated.
- Cleaning and user maintenance must not be performed by children.
- Keep children under 8 years of age away from the toothbrush.
- Do not leave children unattended while using the toothbrush.
- Do not allow children to play with the toothbrush. Children could swallow small parts and choke.

### **Risk of injury**

Handling the toothbrush improperly may result in injury.

- Lay the main cord so that it does not pose a tripping hazard.
- Immediately replace worn or damaged brush heads to prevent injury to the mouth and tooth area.

### **Risk of damage**

Improper handling of the toothbrush may damage it.

- Do not drop the hand piece. If this does happen, do not use the toothbrush any more. Have it checked by an authorised service provider.
- Turn the toothbrush off when you change the brush head.
- Place the charging base on an easily accessible, level and dry surface.
- Ensure that the main cord, the charging base and the toothbrush do not come into contact with hot surfaces or open flames.
- Stop using the toothbrush if its plastic parts exhibit cracks or breaks or are deformed.

### **First Use**

Checking the toothbrush and package contents.

### **Risk of damage**

If you are not cautious when opening the packaging with a sharp knife or other pointed object, you may damage the toothbrush and injure yourself.

- Be very careful when opening the packaging.

1. Take the toothbrush and its accessories out of the packaging.
2. Check to make sure that all the items are included.
3. Check whether the toothbrush or the individual parts are damaged. If this is the case, do not use the toothbrush. Contact Ordo using the customer service contacts provided.

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<b>Problem</b>	<b>Possible Cause</b>	<b>Solution</b>
Lo spazzolino non si carica.	La base di ricarica è difettosa.	Sostituire la base di ricarica entro il periodo di garanzia. Contatta il nostro servizio clienti.
The toothbrush operates with only a small amount of vibration.	The rechargeable battery is depleted.	Charge the toothbrush for at least 8 hours.
The toothbrush does not run for 2 minutes.	The rechargeable battery is depleted.	Charge the toothbrush for at least 8 hours and make sure that the toothbrush is properly placed on the charging base.

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Thanks for Joining

# ORDO

Welcome to the  
Oral Care Revolution

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


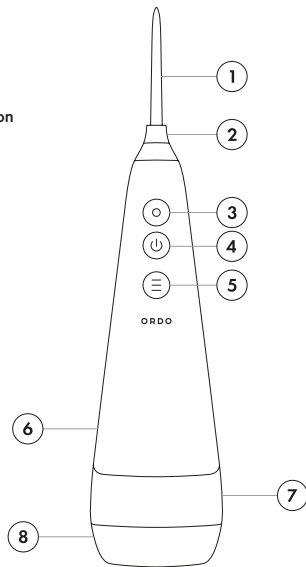
9. Mode indicator panel
10. Floss tip release button
11. USB-C charging port
12. Water tank release button
13. Water tank cap
14. Water tank base release button

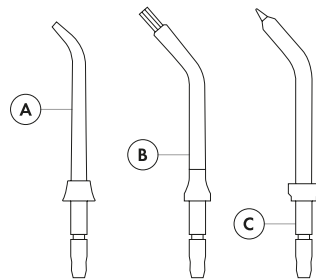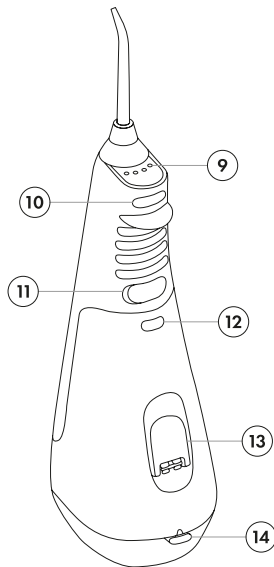
- A. **Premium Clean** – Ideal for excellent everyday cleaning. Suitable for most people.
- B. **Orthodontic** – Tapered bristles help clean around braces, bridges, implants and other dental work.
- C. **Periodontal** – Designed to rinse periodontal pockets, helping to protect against gum disease.


Orthodontic & Periodontal tips sold separately. Our full selection of floss tips are available at [ordolife.com](http://ordolife.com)

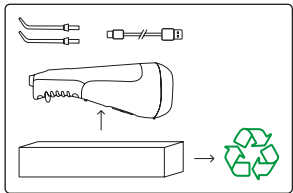**(1)** Remove all the contents from the box and check you have all parts specified. Make sure you recycle the inner and outer boxes after use.

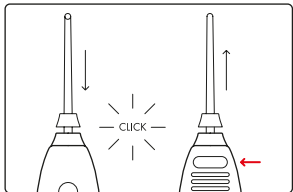
**(2)** Simply push the floss tip into the slot at the top of the water flosser, until it clicks into place. To eject, press the floss tip release button on the back of the water flosser, then pull the floss tip out.

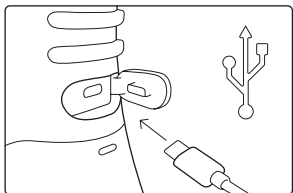
**(3)** Before using your water flosser, you should charge it for 6 hours. Open the USB-C port cover, then plug in the USB-C connector of the charging cable. Connect the USB-A connector to a power source. Ensure the water flosser is not on prior to charging.


**MAKE SURE THE USB-C PORT IS DRY BEFORE CHARGING!**

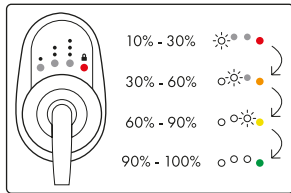
**(4)** During the charging cycle, the lights on the mode indicator panel will follow a sequence as shown in the diagram, to indicate the water flosser's charging status. When the coloured light turns green, the water flosser is fully charged. A full charge will last 30 days when using the water flosser for 2 minutes each day.

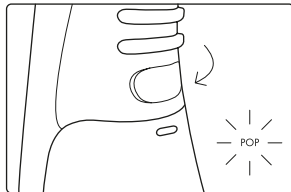
**(5)** Once you have charged the water flosser, remove the USB-C charging cable and securely close the cover. Make sure that this is closed properly to avoid water getting inside.

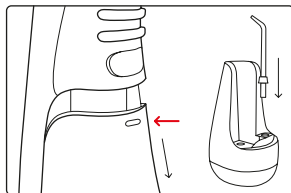**(6)** For travel convenience, you can store your Premium Clean floss tips inside the water tank, as shown in the diagram. Simply push the water tank release button at the top of the water tank and slide off. Inside you will find two slots to securely store your Premium Clean floss tips. To reattach the tank, slide back onto the main body until it clicks into place.

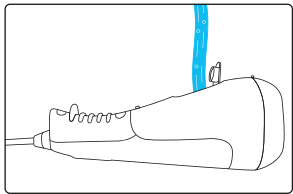
**(1)** To fill your water flosser, flip open the cap on the back of the water tank and fill with clean, lukewarm water (under 40°C). You could also add a few drops of Ordo Mouthwash Concentrate for extra freshness. If you use mouthwash, run a full tank of fresh water through the flosser after use, to avoid residue build-up.

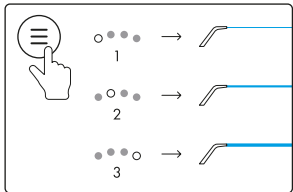
**(2)** Short press the mode selection button to scroll through the 3 pressure settings:

- 1 - Sensitive Mode (low) – recommended for new users.
- 2 - Daily Clean Mode (medium)
- 3 - Power Clean Mode (high)


You can change modes whilst using your water flosser.

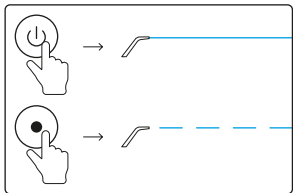
**(3)** Lean over the sink, place the floss tip in your mouth and part your lips slightly to allow water to flow out. Press the power button to begin flossing. Alternatively, press the on-demand button to control the flow of water – press and hold for a continuous flow and release to stop.

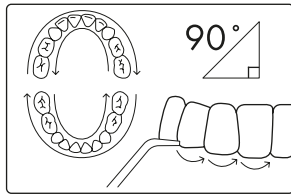
**(4)** Hold the floss tip at a 90° angle to your gum line and pause in between each tooth to ensure a thorough clean. Clean along the gum line as well as between teeth. Clean inside as well as outside surfaces. To get the most from the water tank, do not tilt the water flosser back during use. To help protect your water flosser, it will automatically turn off after running for 2 minutes.

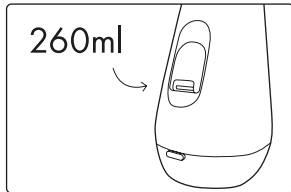
**(5)** Your water flosser has a 260ml tank capacity, which provides 55 seconds of use. We recommend 2 full tanks per flossing session, one for the top teeth and one for the bottom. Make sure the tank is fully emptied after use and open the cap to let it air dry. Wipe down the floss tip and main body.

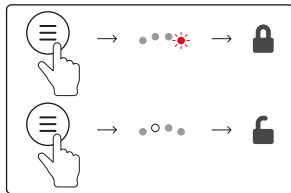
**(6)** To deactivate your water flosser for travel or storage, press and hold the mode selection button until the red light begins to flash, at which point the Lock Mode is activated. In Lock Mode the lock light will flash 5 times when you short press any button. Simply press and hold the mode selection button again to unlock the device – the last used mode will illuminate once your flosser is unlocked.

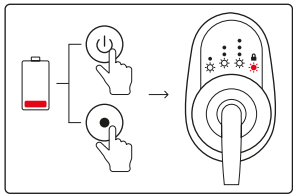
**(1)** At 30% battery, all the lights on the mode indicator panel will flash 5 times after use. When the water flosser has no power, the lights will also flash 5 times, but the device will not turn on – you need to recharge to a sufficient level before you can continue using the water flosser.

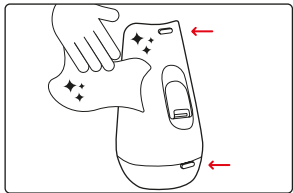
**(2)** To clean the water tank, first remove it from the main body by pressing the release button at the top of the tank and sliding off. Remove the base from the tank by pushing the release button at the bottom and pulling open. You can clean all parts of the tank with warm, soapy water or place in the dishwasher.

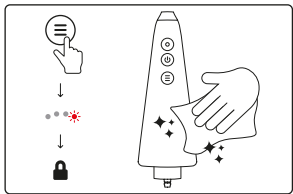**(3)** To clean the main body of the water flosser, ensure it is switched off and disconnected from any power source. We recommended activating the Lock Mode for cleaning. To avoid scratching, use a soft, dry or damp cloth to clean the main body.

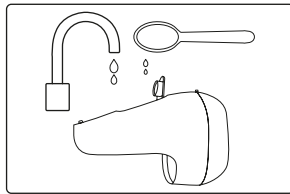**(4)** To clear hard water deposits from within your water flosser, fill the tank with water and add 1 tablespoon of white vinegar. Run the water flosser until the tank is empty, then rinse thoroughly with clean, warm water. Repeat this process every 1 to 3 months, depending on the hardness of your water.

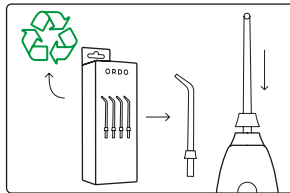**(5)** It is recommended that you replace floss tips when they become worn or dirty, this is usually every 6 months. Replacement tips as well as Orthodontic and Periodontal tips are available at [ordolife.com](http://ordolife.com)

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