CF-311+

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Directory

Product Introduction	09
Functions	10
Applications	11
Specifications / Operation Guidelines / Product Maintenance	12
Attention / Fault Diagnosis	13

MEDIVON INFRARED BODY SLIMMING MASSAGER

Product Introduction

Infrared Body Slimming Massager is exclusively developed by a professional design team, which based the design on human physiology and combined those ideas with modern work and life habits. It has a simple appearance and is easy to carry and operate. Four different massage heads with different functions on different parts. The high frequency vibrations smash cellulite on multiple sites. Infrared hot compress accelerate fat burning. You could enjoy it anytime and anywhere.

Do not use the massager without the safety shield protection except while using the Ball Head and Frosted Head.

The safety shield prevents hair from being tangled in the head. It is unsuitable for ball-type and frosted-type massage heads.

MEDIVON INFRARED BODY SLIMMING MASSAGER

Functions

Fat Buster — Infrared Body Slimming Massager

Vibration Massage to Relieve Fatigue

Powerful frequency vibration 3500 times/minute dredge meridians. Targets directly into pain acupoints to relieve deep muscle strain.

Portable and Powerful Body Shaping

This is a portable and powerful body shaping product designed to promote blood circulation, accelerate fat burning, and improve body immunity. It is a healthy way to lose weight safely.

Improve Complexion, Promote Blood Circulation

High frequency vibrations stimulate and massage body parts, promote blood circulation, and help you enjoy a healthy lifestyle easily.

Deep Detoxification, Beautify Bottom Line

High frequency vibration helps to clear stool and improve constipation for detoxification and positioning massage to shape a sexy bottom line.

MEDIVON INFRARED BODY SLIMMING MASSAGER

Applications

Neck and Shoulder


Arm


Abdomen

Waist

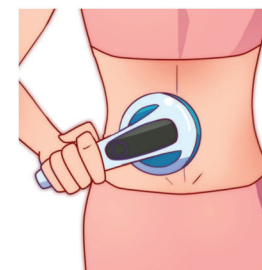
Back

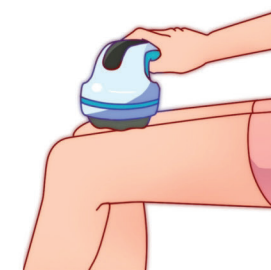Thigh

Calf

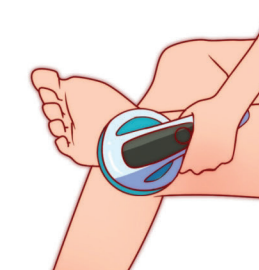Foot

EASY · COMFORTABLE · DEXTEROUS · SIMPLE

MEDIVON INFRARED BODY SLIMMING MASSAGER

Specifications

Product Name	Infrared Body Slimming Massager	Net Weight	0.9 Kg
Product Model	CF-311+	Size	20*17*10 cm
Voltage	AC 220-240V	Materials	ABS/PP/TPV
Power	25W	Standard	GB4706.1-2005
Gross Weight	1 Kg		GB4706.10-2008

Operation Guidelines

- 1.Insert the plug into a power socket.
- 2.Rotate power dial to the right to start the machine.
- 3.Adjust the intensity levels by rotating the dial to the right or left.
- 4.It is recommended to use lower massage strength for initial use and then gradually increase intensity.
- 5.Turn the device off by rotating the button to utmost left, then unplug.
- 6.If initial use, massage through the clothes. Use direct contact with the skin after several times of use.

Product Maintenance

- 1.Use soft cloth o cleanser to scrub the machine. Removable massage heads could be cleaned with water individually.
- 2.Do not submerge this machine in water or any other liquid.
- 3.After use, store in a dry place avoiding high temperature, direct sunlight and corrosive gas.
- 4.After use, make sure pull out the plug from socket. Do not pull the power cord to unplug the device.

MEDIVON INFRARED BODY SLIMMING MASSAGER

Attention

- 1.Use within one hour after a meal is not recommended.
- 2.Pregnancy and menstruating women, skin disease patients, persons with wounds, and heart disease patients should not use this product.
- 3.Children or persons with disabilities should be supervised.
- 4.Do not use product in the bathroom, swimming pool, or other damp or wet places.
- 5.If discomfort is experienced, stop using the product immediately.
- 6.Do not massage in the same place for more than 15 minutes.
- 7.Reduce usage times, as it can prolong service life.

Fault Diagnosis

Fault	Ireason	Troubleshooting
The motor do not work?	The plug is not inserted firmly.	Insert the plug firmly.
	Switch is not in place	Check the switch.
Not Work after a long time working.	Thermal protection	30 minutes rest

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Medivon Sp. z o. o.
ul. Solec 18, lokal B-21,
00-410 Warszawa, Poland

English

Hereby, Medivon Sp. z o. o. declares that the device equipment Medivon CF-311+ is in compliance with Directive 2014/53/EU. The full text of the EU declaration of conformity is available at the following internet address:

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